


A Roadmap to Successful Flash Chromatography

July 26, 2018

Bob Bickler



- » Introduction to flash chromatography
- » Where to begin?
- » Normal-phase flash chromatography
 - » Method development and optimization
 - » Sample loading options
- » Reversed-phase flash chromatography
 - » Method development and optimization
 - » Lipophilic compounds
 - » Ionic/ionizable compounds
- » Scaling-up

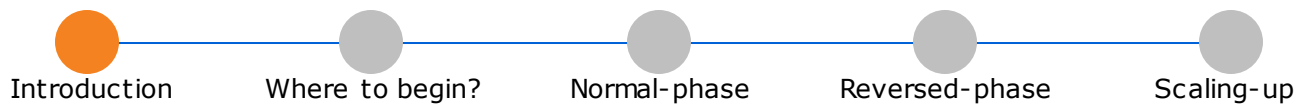

What is flash chromatography?

- » Flash chromatography is a preparative liquid chromatography purification technology that separates and isolates one or more compounds from a mixture of compounds

How does it work?

- » It works by using...
 - » A solid sorbent or media (stationary phase) to attract the sample's compounds
 - » A solvent or solvent blend (mobile phase) to *selectively* elute the compounds from the stationary phase

Basic terminology


- » **TLC** = Thin-layer chromatography. Used to evaluate solvent mixtures for flash method optimization
- » **R_f** = Retardation factor. A measure of a compound's elution on a TLC plate relative to the solvent migration. Used to help determine how much solvent is required to elute the compound (expressed in CV)
 - » $R_f = D_1 / D_{SF}$
 - » Scale is 0-1
- » **CV** = Column volume. The volume of a packed flash column or cartridge specified in mL. Also used as a measure of compound retention, $CV = 1/R_f$
- » **ΔCV** = Delta CV. The difference in elution volumes for two adjacent compounds. This value is used to determine loading capacity for a mixture of compounds
- » **Selectivity**. The amount of separation between adjacent eluting compounds
- » **Isocratic**. A mobile phase in which the solvent ratio remains constant throughout the purification
- » **Gradient**. A mobile phase that changes polarity (a.k.a. strength) over the course of the purification moving from a weaker solvent blend to a stronger solvent blend. There are two gradient types commonly used
 - » Linear gradient. A gradient that is created at a constant rate (% strong solvent/CV)
 - » Step gradient. A gradient created by increasing solvent strength in a step-wise fashion.
- » **Equilibration**. The process used to prepare the flash column for purification. Equilibration displaces the air in a new column and conditions the stationary phase to a constant polarity. Typically, 3 to 5 CV are required for full equilibration
- » **Dry load**. A sample introduction technique in which the sample is dried onto a sorbent prior to putting inline with the flash column. The typical benefit is improved separations and better fraction purity
- » **Flow rate**. The solvent pumping rate in mL/min
- » **Linear velocity**. Solvent flow expressed in distance per unit of time (mm/sec or cm/min). Used when scaling up a purification to ensure a consistent sample/sorbent interaction rate (kinetics)

Flash Chromatography

Elution Methods


» Isocratic

- » Two solvents are mixed at a specific ratio (based on TLC) and kept at that ratio during elution
- » Compound elution typically matches TLC prediction based on the relationship $CV=1/R_f$

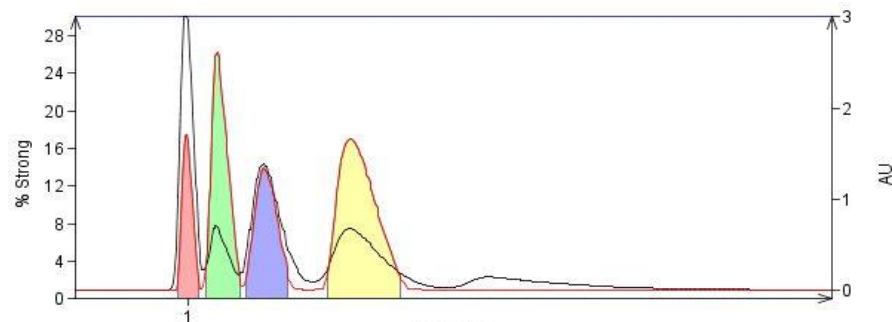


» Linear gradient

- » Two independent, miscible solvents are pumped and blended at continually changing ratios prior to reaching the cartridge
- » Provides better peak shape, improved resolution, improved loading capacity, and higher concentration fractions

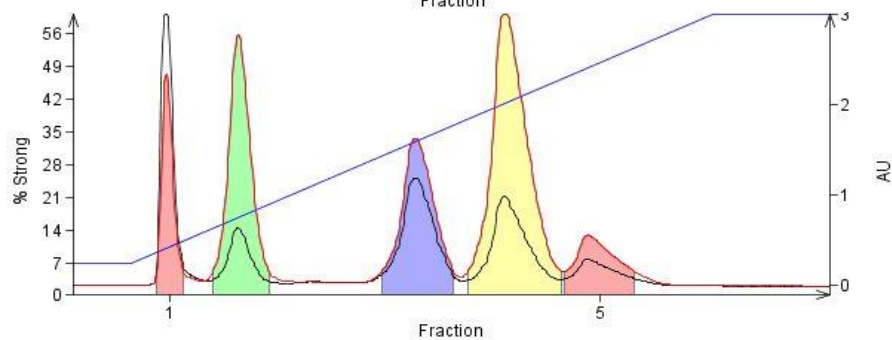

» Step gradient

- » Two independent, miscible solvents are pumped and blended step-wise at increasingly stronger ratios prior to reaching the cartridge
- » Can provide faster elution, increased sample loading, enhanced selectivity, improved resolution, higher concentration fractions, and reduced solvent consumption

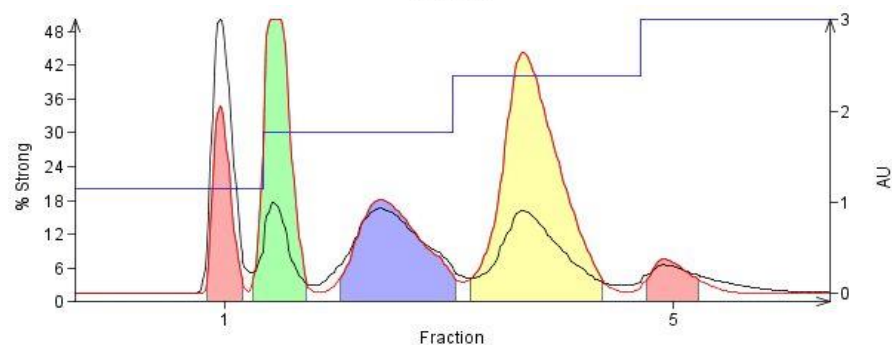

Normal-phase

- » Polar stationary phase
 - » Silica
 - » Polar bonded
 - NH
 - Diol
 - CN
- » Non-polar solvents
 - » Hexane/heptane
 - » Ethyl acetate/acetone
 - » DCM
 - » Alcohols
 - » Ether
 - » Toluene

Reversed-phase


- » Non-polar stationary phase
 - » C18
 - » C4
- » Polar solvents
 - » Water
 - » Methanol/ethanol
 - » Acetonitrile
 - » Acetone
 - » Tetrahydrofuran (THF)

Remember, for normal- and reversed-phase...

Like Likes Like

Where to Begin?

Understand your Purification Goals

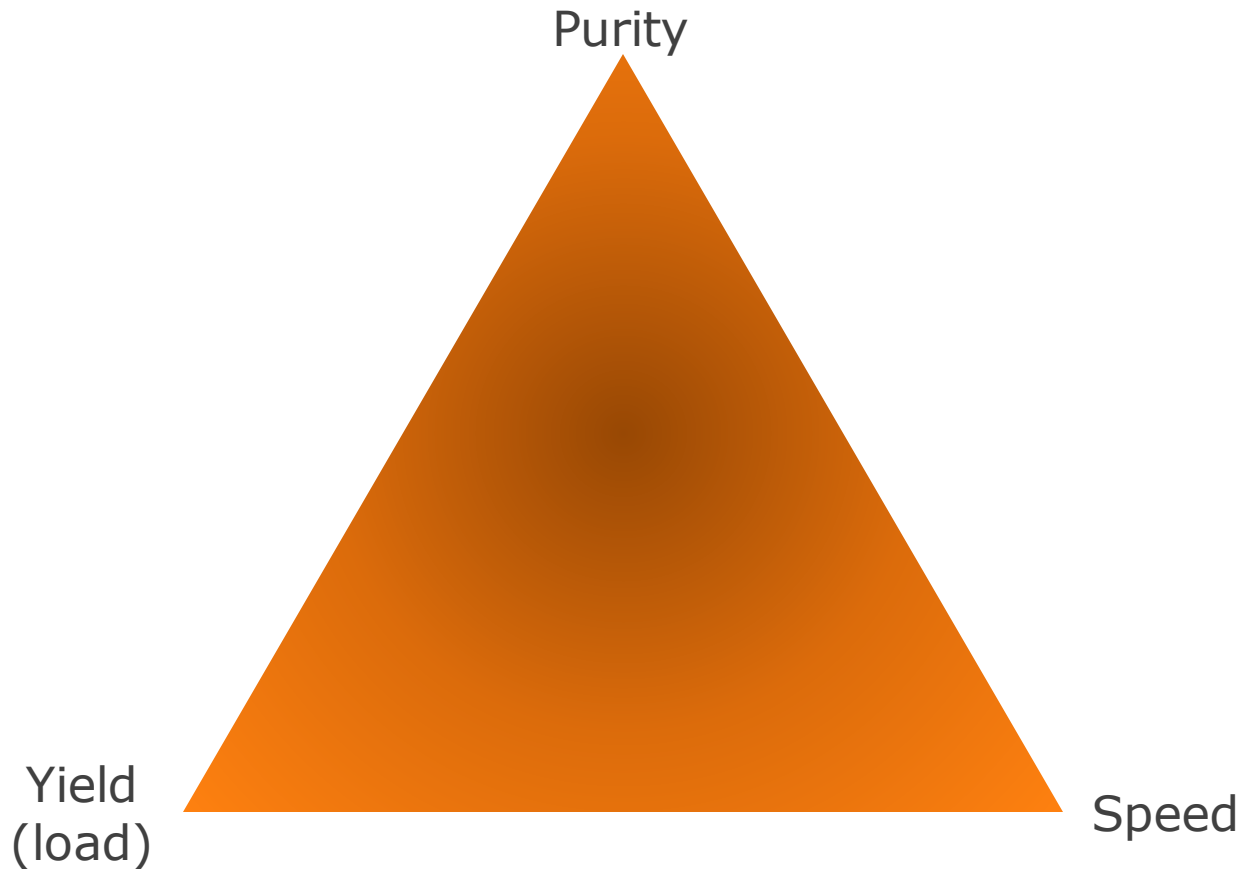

Where to Begin?

Understanding Your Sample

- » Only C and H
 - » Organic soluble compounds
 - » Reversed-phase
 - » Normal-phase

- » Organic amines
 - » Organic soluble
 - Amine-bonded silica
 - » Water soluble
 - Reversed-phase

- » Most others
 - » Normal-phase silica (if organic soluble)
 - » Reversed-phase
 - Water soluble or soluble in higher polarity solvents
 - Ionic/ionizable
 - Acidic/basic

Where to Begin?


Understanding Your Sample - Solubility

Rank	Solvent polarity	Mode
1	Dimethyl sulfoxide (DMSO)	Reversed-phase
2	Dimethyl formamide (DMF)	Reversed-phase
3	Water	Reversed-phase
4	Methanol (MeOH)	Reversed-phase
5	Acetonitrile (MeCN)	Reversed-phase
6	Ethanol (EtOH)	Reversed-phase
7	Propanol (PrOH)	Normal-phase*/Reversed-phase*
8	Acetone	Normal-phase*/Reversed-phase*
9	Tetrahydrofuran (THF)	Normal-phase*/Reversed-phase*
10	Ethyl acetate (EtOAc)	Normal-phase*/Reversed-phase*
11	Ether	Normal-phase/Reversed-phase
12	Methylene chloride (DCM)	Normal-phase/Reversed-phase*
13	Toluene	Normal-phase/Reversed-phase*
14	Hexane	Normal-phase/Reversed-phase*

* Dry loading is best for these solvents

Method Development

Normal-phase

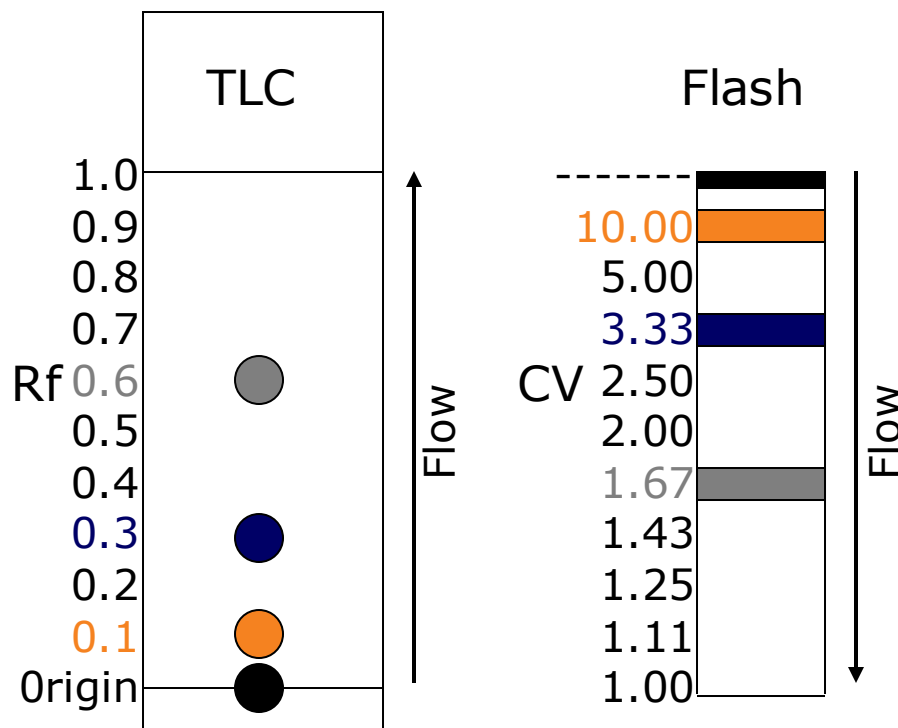

- » Compound retardation factor (R_f)
 - » Inversely proportional to flash elution volume (CV) in isocratic elution mode
 - $CV = 1/R_f$
 - » ΔCV used to determine sample load

Method Development

Normal-phase Elution Solvents

Solvent	Group	Strength
Methanol	II	0.70
Ethanol	II	0.65
Propanol	II	0.60
Tetrahydrofuran	III	0.53
Acetone	VIa	0.50
Acetonitrile	VIb	0.50
Ethyl acetate	VIa	0.43
Ether	I	0.40
Dichloromethane	V	0.32
Toluene	VII	0.22
Heptane	---	0.01
Hexane	---	0.01

- » Scout solvents from different selectivity groups
 - » Maximize separation (ΔCV)
- » Adjust solvent strength to elute target within the R_f range of 0.15 to 0.4


Method Development

Normal-phase



	<u>Rf</u>	<u>ΔRf</u>	<u>CV</u>	<u>ΔCV</u>
	0.10		10.0	
Optimal Range	0.15		6.7	
	0.20	0.05	5.0	1.7
	0.25	0.05	4.0	1.0
	0.30	0.05	3.3	0.7
	0.35	0.05	2.8	0.5
	0.40	0.05	2.5	0.3
	0.45			2.2
0.50			2.0	
0.60			1.6	
0.70			1.4	
0.80			1.25	
0.90			1.11	
1.00			1.0	

R_f – CV relationship



ΔR_f or ΔCV?

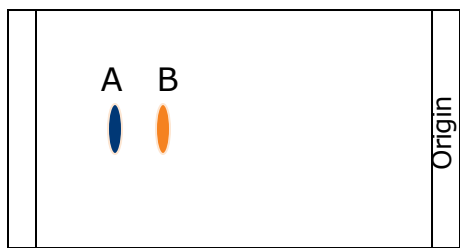


R_fA = .80 R_fB = .67 ΔR_f = .13
ΔCV = 0.24

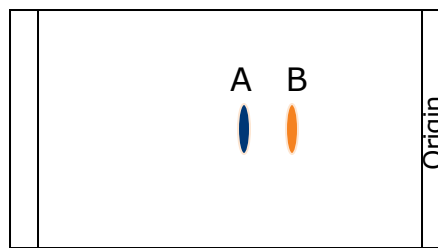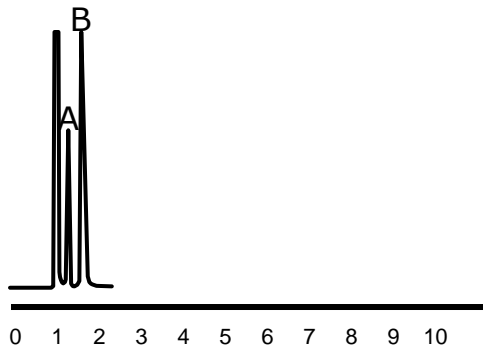

R_fA = .47 R_fB = .34 ΔR_f = .13
ΔCV = 0.81

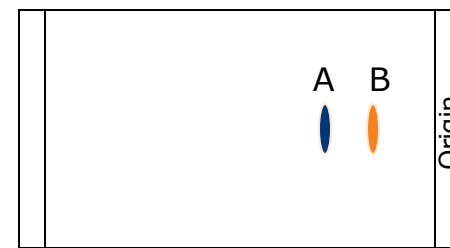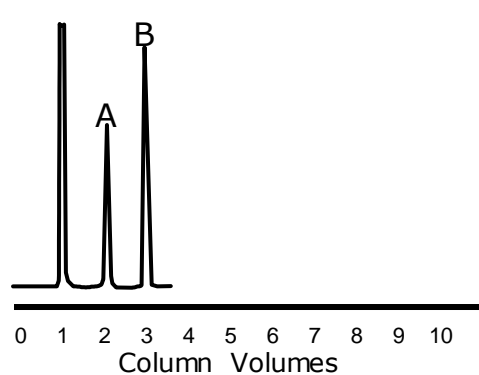


R_fA = .32 R_fB = .18 ΔR_f = .14
ΔCV = 2.43

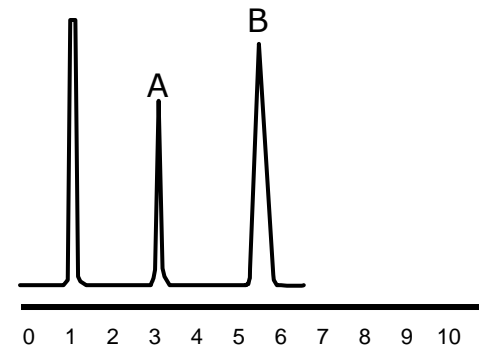

- » ΔR_f constant
- » Increasing ΔCV with decreasing R_f
- » Predict maximum sample loading with ΔCV, not ΔR_f


Solvent Choice Impact

Solvent Chemistry Matters

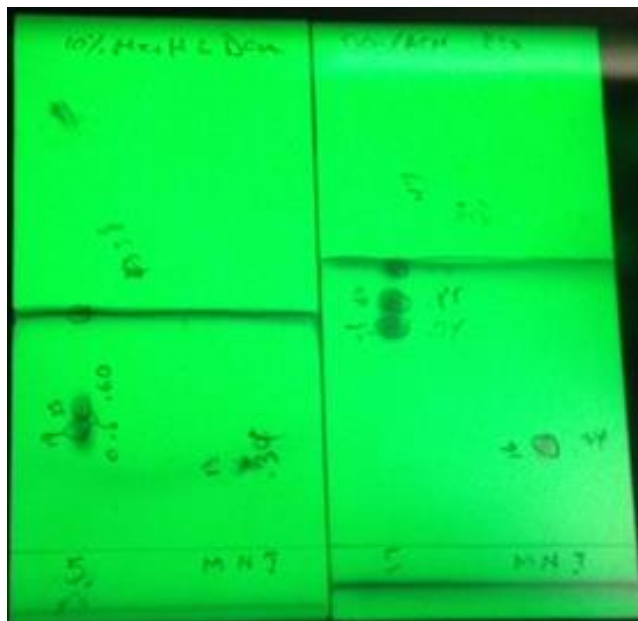

Goals and Chemistry

- » Goal
 - » Separate and isolate each compound
- » Solubility
 - » Alcohol
 - » Acetone
- » Chemistry
 - » Contains carbon, hydrogen, and oxygen
 - Hydroxyl groups
 - Carbonyls

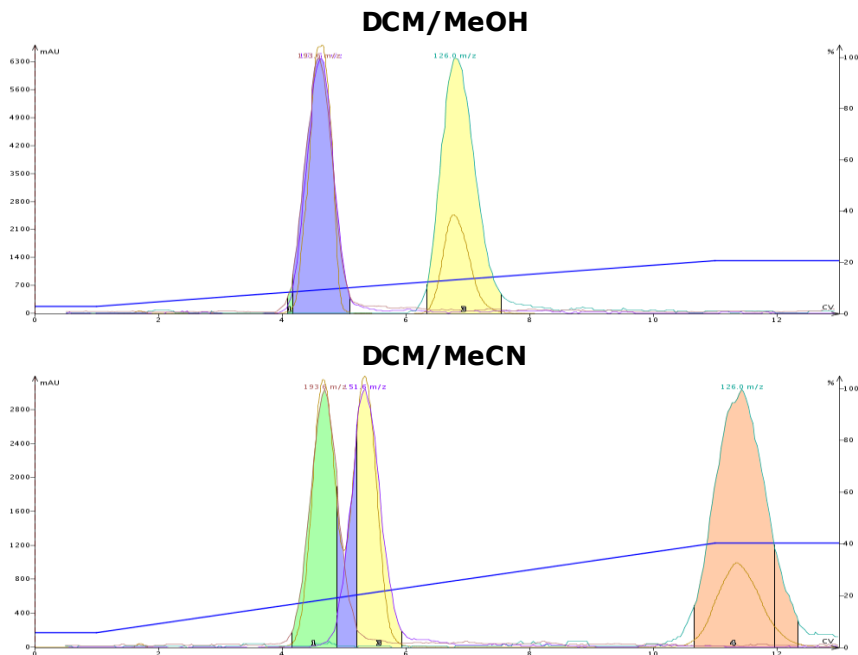


Method Development and Chromatography

- » TLC
 - » DCM/MeOH (9:1)
 - » DCM/MeCN (8:2)
 - » Same solvent strength (0.36)
- » Flash
 - » Liquid load in acetone (150 mg/0.15 mL)
 - » DCM/MeOH unable to separate two compounds
 - » DCM/ACN successful


Loading Techniques

Liquid or Dry?

1-mL liquid load (acetone), 200 mg sample,
10 gram Biotage® SNAP Ultra, hexane/EtOAc gradient

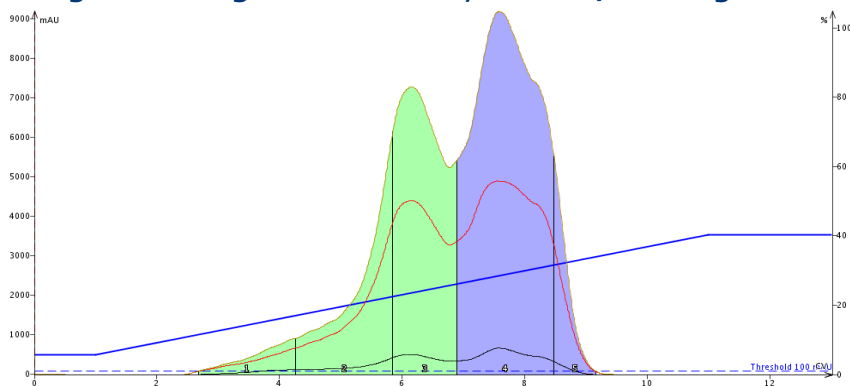

Dry load (1 gram silica), 200 mg sample,
10 gram Biotage® SNAP Ultra, hexane/EtOAc gradient

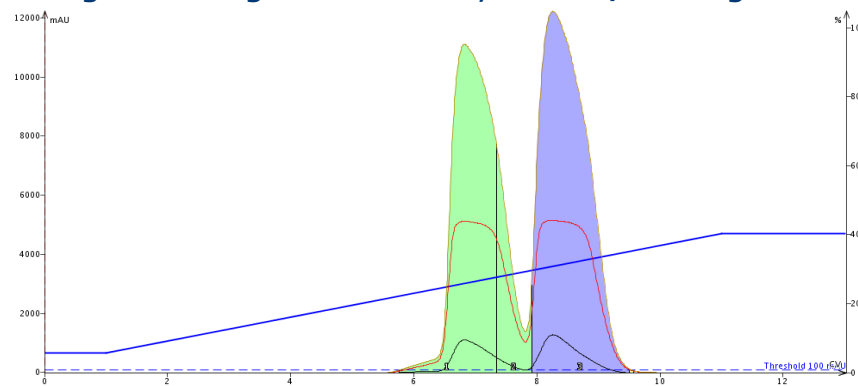

Liquid Loading

- » Weak dissolution solvent
- » High sample concentration
 - » >1 g/mL
- » Low load volume
 - » <1% of a CV


Dry Loading

- » Solvent dissolves sample completely
- » Solvent is volatile
- » Mix with sorbent and dry completely

Dry Loading


Options and Techniques

- » Select a sorbent
 - » Silica
 - » Diatomaceous earth (ISOLUTE® HM-N, Celite)
 - » Alumina
 - » Florisil®
 - » Ion exchanger (silica-based)
 - » C18
- » Weigh enough sorbent
 - » Good sample to sorbent ratios range between 1:2 to 1:4
- » Dissolve sample in minimal volume of volatile solvent and mix with sorbent
- » Dry completely
 - » Free-flowing powder
- » Use empty Samplet® cartridge or dry load vessel (DLV)

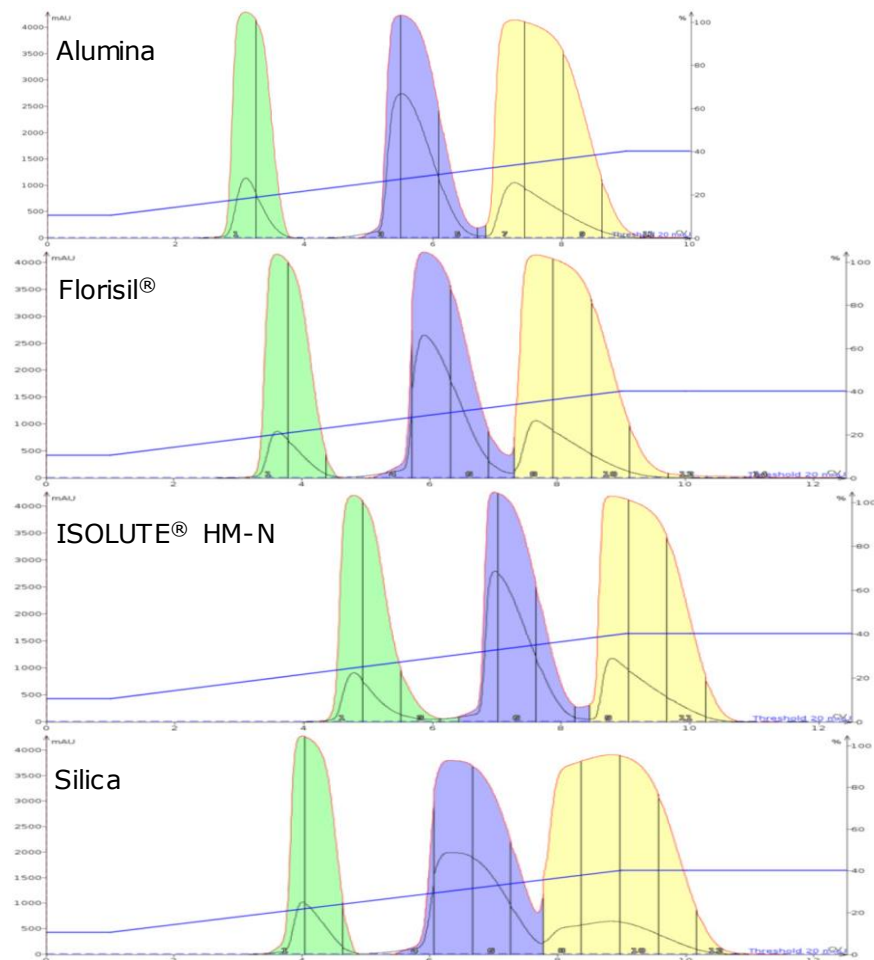

Dry Loading

Options and Techniques

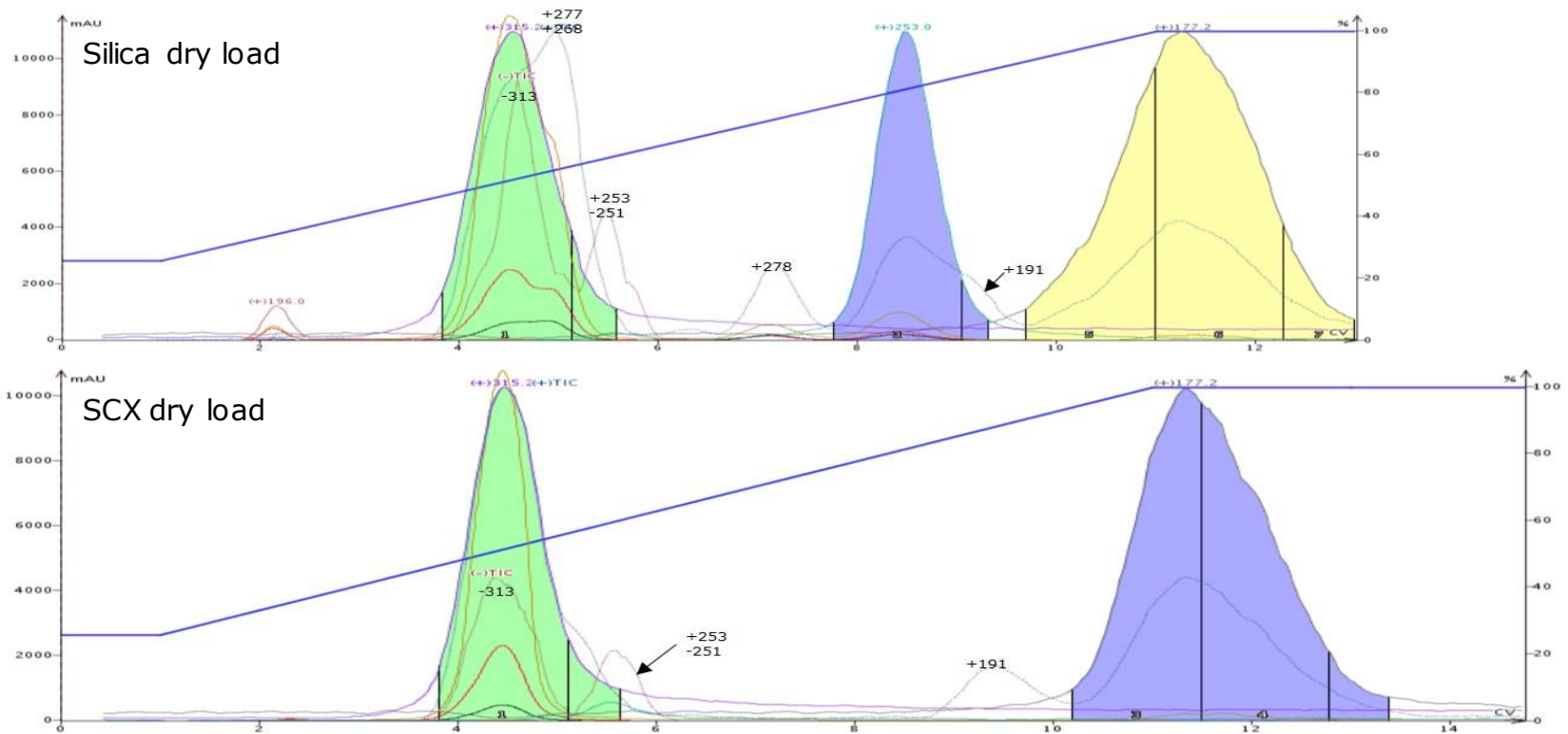


» Impurity scavenging

- » Remove interfering compounds using a scavenger such as an ion exchange media

Flash Chromatography in Practice



Natural Product Purification - Cannabis Extract

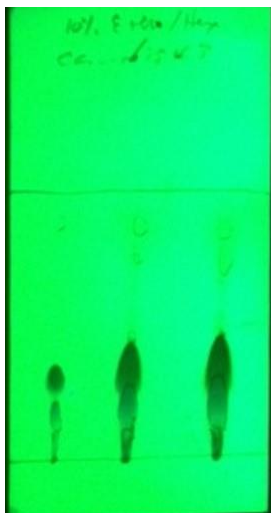

Waxes, terpenes, carotenoids

Cannabinoids

Cannabinoid acids
Chlorophyll and other polars

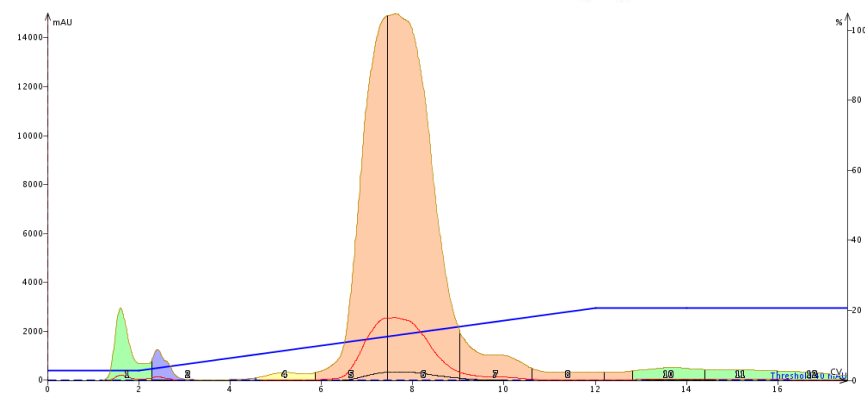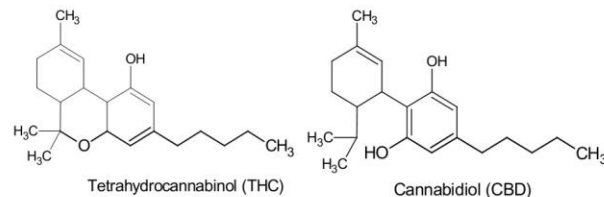


Goals and Sample Chemistry

- » Goals
 - » Clean-up extract to isolate THC + related compounds from undesirable compounds
- » Chemistry
 - » Contains carbon, hydrogen, oxygen
 - » Hydroxyl groups
 - » Solubility (hexane)

Method Development and Chromatography

- » TLC
 - » 10% ether in hexane
- » Flash
 - » Biotage® SNAP Ultra, 10 g
 - » Liquid load
 - 250 mg (0.1-mL hexane)

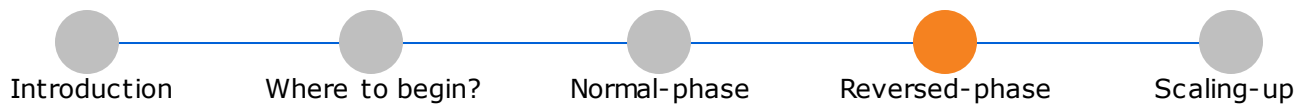

» Use HPLC

» Use a scaling column

- HPLC column packed with reversed-phase flash sorbent

» Adjust solvents and gradient to optimize separation

- Water/methanol
- Water/acetonitrile
- Acetone and THF can be used in place of methanol and acetonitrile for more hydrophobic compounds


» Add volatile acid, base, or buffer as needed

- Improve peak shape
- Improve separation

- » Sample solubility determines starting conditions
 - » Water, DMSO, DMF soluble
 - Start gradient at 5% B for 1 CV
 - B is usually methanol or acetonitrile
 - Finish gradient at 50% B (over 10 CV)
 - » Methanol, acetonitrile, or other organic solvent soluble
 - Start gradient a 50% B for 1 CV
 - Finish gradient at 100% B (over 10 CV)
- » If compound of interest elutes too early, reduce final % B
- » If target elutes too late, increase starting % B
- » If target does not separate from impurities change solvent B, modify pH/ionic strength, or change gradient slope

Flash Chromatography

Cannabis Extract

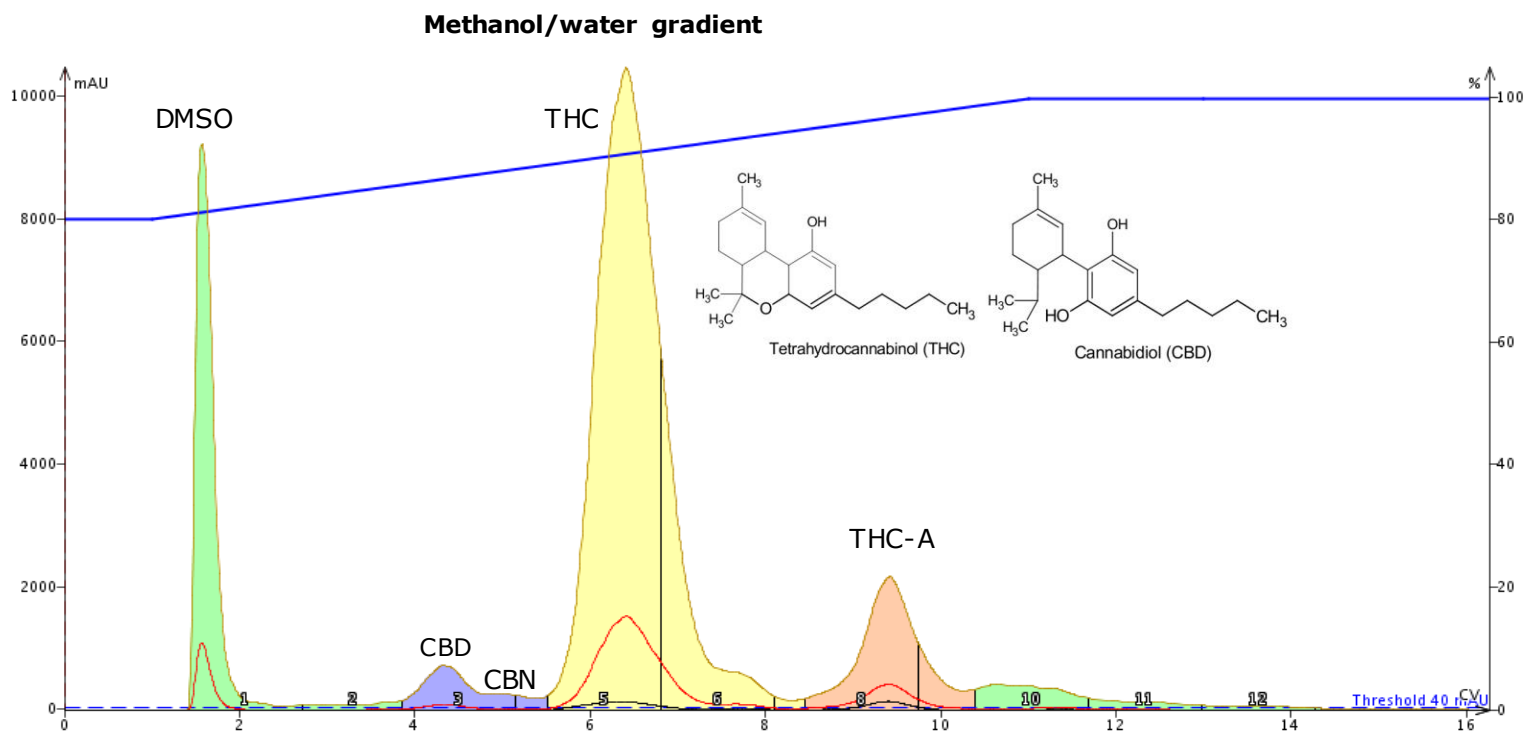

- » Goal
 - » Isolate CBD, CBN, THC, other components


- » Solubility
 - » Hexane
 - » Alcohols
 - » DMSO

- » Chemistry
 - » Contains carbon, hydrogen, oxygen
 - Hydroxyl groups
 - Ether


- » Flash
 - » Biotage® SNAP Ultra C18, 12 gram
 - » High % organic solvent
 - » 50 mg load, DMSO

Flash Chromatography

Ionic, Water-soluble Compounds

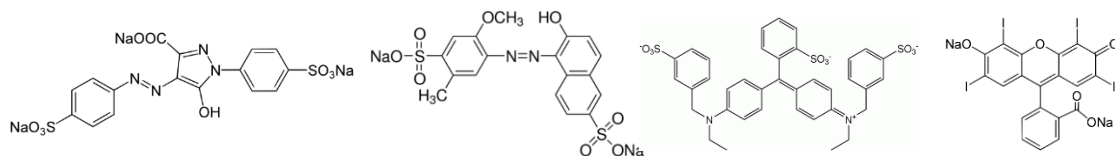


Methanol/water gradient

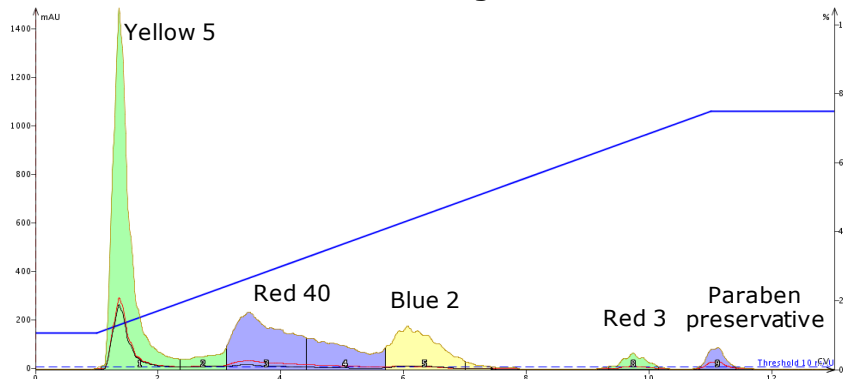


- » Goal
 - » Separate and isolate each dye
- » Solubility
 - » Water
 - » Alcohol
- » Chemistry
 - » Ionic/ionizable
 - » Carboxylic acid groups
 - » Sulfonic acid groups
 - » Quaternary amine

Methanol/water (with 0.1M ammonium acetate, pH 7) gradient

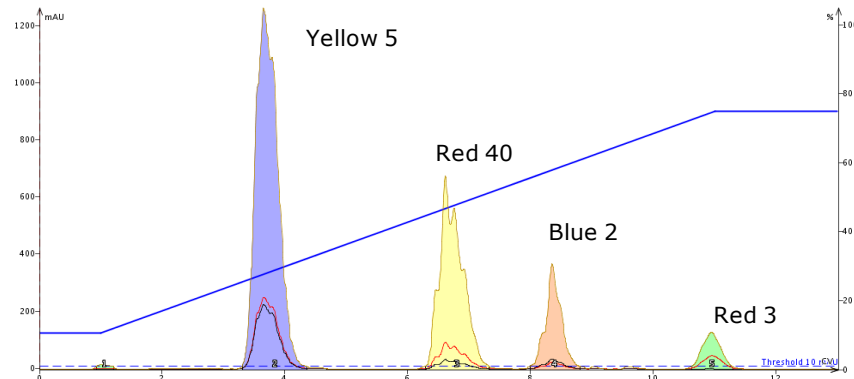

- » Flash
 - » Biotage® SNAP Ultra C18, 12 gram
 - » Low organic solvent
 - » Un-buffered mobile phase
 - Low retention
 - Broad peaks
 - » Buffered mobile phase
 - Increased retention
 - Sharp peaks

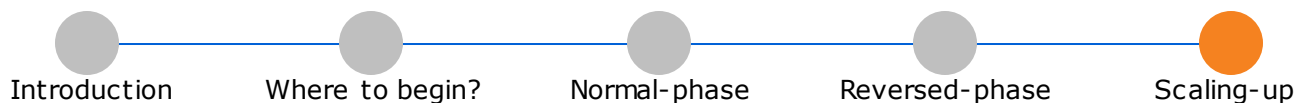


Milligrams to grams to kilograms

- » Small-scale flash easily scaled
 - » Use scale-up factor (SF)
 - » Multiply SF times small-scale load to get large scale load
- » Scale factor equal to the column volume or media mass ratio
 - » Media must be identical
 - » $SF = CV_L / CV_S$
 - » $SF = \text{Large column grams} / \text{Small column grams}$


Flash Chromatography

Scaling up

SNAP KP-Sil 100 g
2.2g load

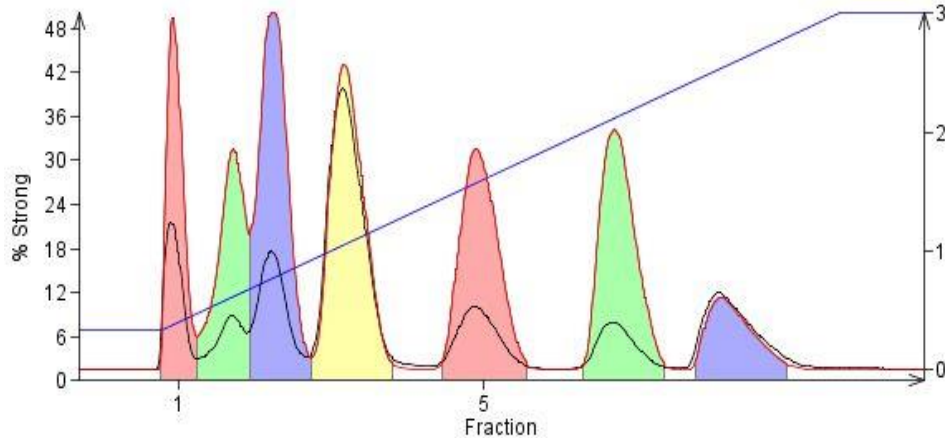

SNAP KP-Sil 340 g
SF=3.4
7.7g load

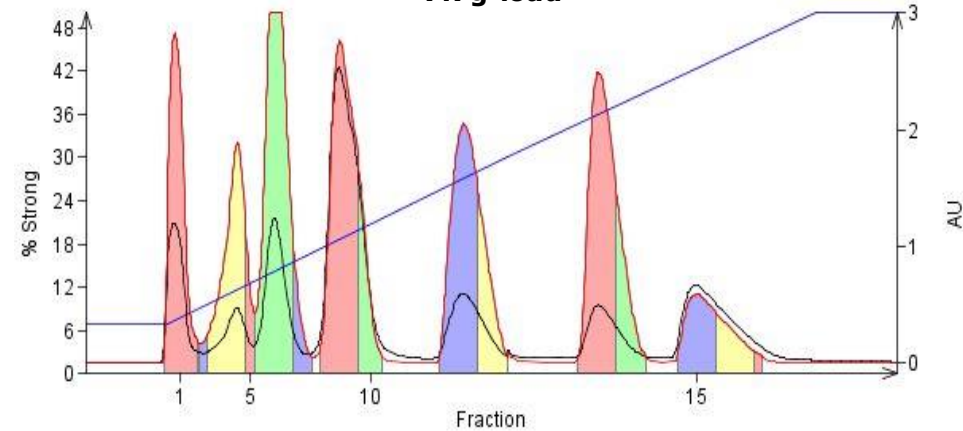

- » If a 10 gram cartridge purified 100 mg, then a 100 gram cartridge will purify a gram and a 1000 gram cartridge will purify 10 grams
 - » Stationary phase must be identical
- » No method modification required other than flow rate and rack size
 - » Flow rate does not scale!
 - » Need to match linear velocity

Flash Chromatography Success

Summary

- » Flash chromatography
 - » Determine purification goal(s)
 - » Determine crude sample solubility
 - Liquid or dry load
 - » Understand possible functional group impact on the purification
 - Normal- or reversed-phase
 - Neutral or modified solvents
 - » Spend a little time on method development